


ADPT-NN-OM32-STIM

Stimulation Adapter for MPA32I and NeuroNexus Probe with Omnetics Connector

Reference REF


The reference is directly connected to ground. Other connections are possible on request.

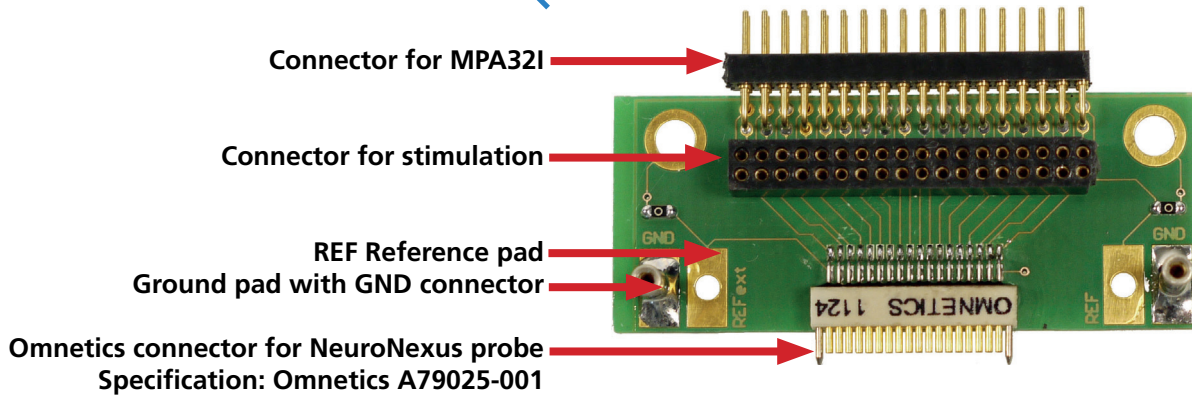

Please connect the MPA32I with the bottom side up to the top side of the ADPT-NN-OM32.

- Pin 1: Ground
- Pin 2: Reference
- Pin 3 to 34: Channel 1 to 32
- Pin 35: Reference
- Pin 36: Ground

Pin Layout of the ADPT-NN-OM32 (Top Side)

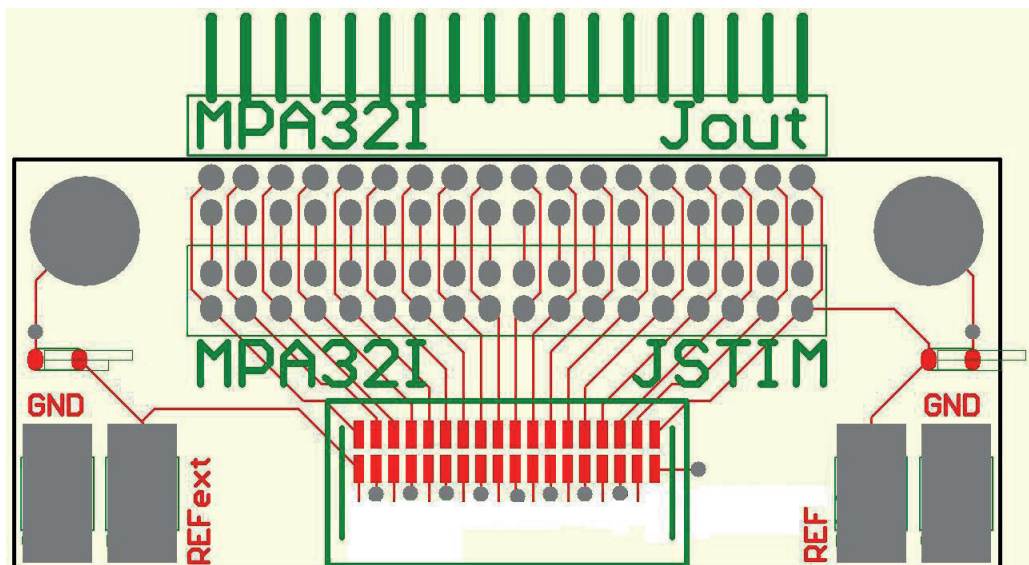

June 2021