

ADPT-W2100-HS8/16-90° for Use with the W2100-System



Adapter to connect a W2100-HS8/16 Headstage rectangular to the Surface

Applications


Using two of these rectangular adapters ADPT-W2100-HS8/16-90° gives the possibility of using two W2100-HS8/16 headstages at once.

Connect both headstages via adapters vis à vis at each other and record from 16 or 8 electrode channels simultaneously.

ADPT-W2100-HS8/16-90°
Pin layout of the upright Omnetics Connector
Orientate the headstage as shown on the picture.
Please read also the datasheets of the W2100-System headstages W2100-HS16 or W2100-HS8.

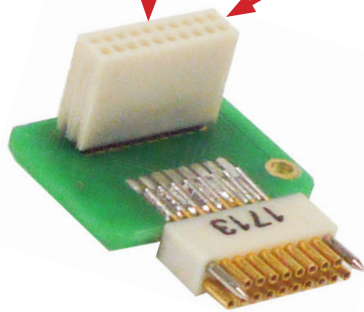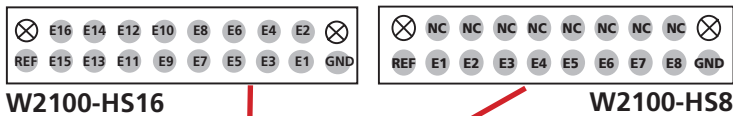

ADPT-W2100-HS8/16-90°

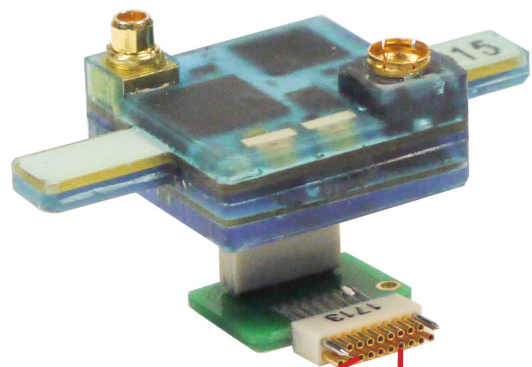

**ADPT-W2100-HS8/16-90°
connected to a W2100-HS16**


W2100-HS16

⊗	Guide post
E1 to E16	Recording electrodes
GND	Ground
REF	Reference


W2100-HS8

⊗	Guide post
E1 to E8	Recording electrodes
GND	Ground
REF	Reference

W2100-HS16

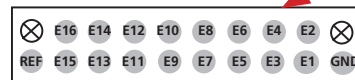

W2100-HS8

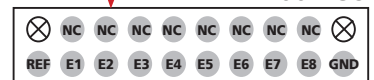

ADPT-W2100-HS8/16-90°

Pin layout for the rectangular Omnetics Connector
Please orientate the headstage as shown in the picture.