


## W2100-System Power Bank: W2100-PB

### Cable Connection to power a Wireless Headstage with a Powerbank for Use with the W2100-System

#### Technical Specifications

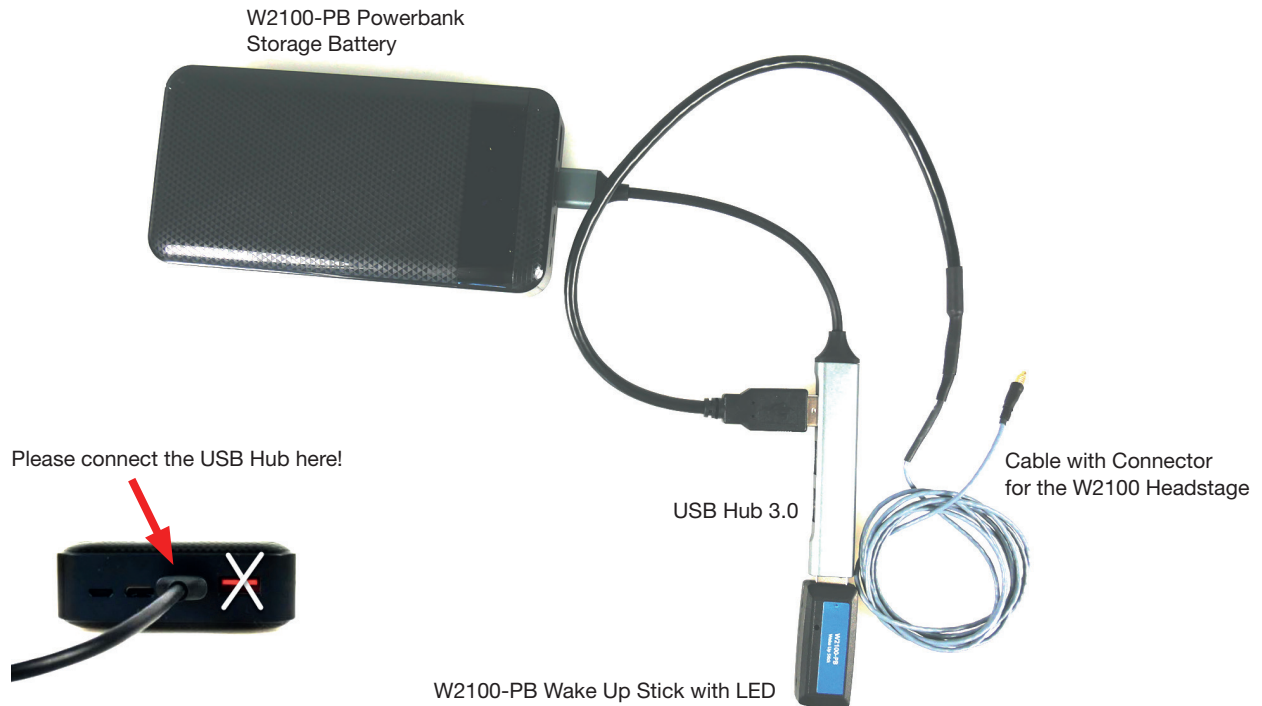

Power Bank	
Powerbank Storage Battery	The powerbank provides (for example 26.000 mAh) to serve a W2100 headstage in a long time experiment without changing the battery.
Charging Set	Charging Set for the Powerbank
USB Hub	Please connect the USB Hub to the USB 3.0 port in the center of the powerbank!
Application	
Power Supply Cable with Headstage Connector	Power supply cable connection to the W2100 headstage /W2100-USB-A-MMCX).
W2100-PB Wake Up Stick	<p>Delivers “Wake up” impulses for the powerbank storage battery. The power consumption of a headstage is very low, that is why the storage battery stops working after a while and the headstage would not be powered any longer.</p> <p>The W2100-PB wake up stick consumes a small amount of power once every 10 secondes indicated by the red LED. These impulses prevent that the powerbank storage battery is switching off.</p>
Note:	Please do not remove the W2100-PB Wake Up Stick from he powerbank during the experiment!